



max steam®

by max stir since 1966

FINISHING



CHECK
ONE BY ONE

QUALITY

DUPLICATION



SERVICE

MODELLING




SPIRIT



EXPERIENCE

MODEL PREPARATION

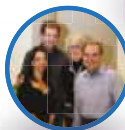

PARTNERSHIP



MASTER SPECIALIST

WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966



Are you looking for an expert and specialized manufacturer?



Max steam is the answer: an Italian family company which develops and makes quality lab equipments since 1966.

Our steam cleaner leads the market as quality and performance, available in manual or automatic way and with different capacity tank. We can offer a wide range of Lab. Machines and innovative equipments from

PREPARATION MODEL → MODELLING
DUPLICATION (complete dentures and implant, fix and removable prosthetics) → THERMOFORMING
→ SANIFICATION & STERILIZATION

Family trade, combination of particular care for details and quality with innovative and flexible vocation to improve the products. Max steam spirit is based on experience, service spares parts and selling support. So the main goal is the relationship with our clients: dealers, dentists and technicians.

Our satisfaction is in your satisfaction!!!

The Maximum quality and service value.

All products are tested and checked individually by our staff and CE conformed. We are involved to energy saving policy.

**WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER**

max steam
LAB MANUFACTURER since 1966



Family trading company and with same spirit family products range.


Every day... hand and mind are cooperating and develop the business together machines, equipments and service solutions.

Our sprit is to listen and join our clients, create a simple, reliable and universal equipments



MS STEAMER had an electronic system to set your favourite pressure

MS MAXFLEXINJECTION also have

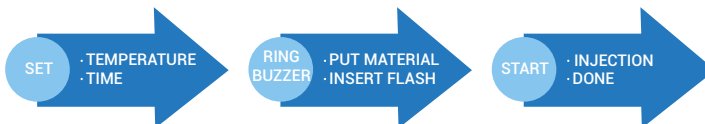

Every system to work with right and satisfaction **SPiRiT**

A brand means much more than...



M

MAXIMUM MAX vision, born with us since 1966. All to be a maximum

MAKER we are real manufacturer

MOVE UP improving and expanding our range

S

SERVICE to deliver spares, to repair, every day

SPECIALIST since 1966

SYSTEM wide range of italian equipments, a full system of quality device for the lab

Every products is check one by one. Made in Italy. All Italian components. Every model is available at 230V or 110V. The equipments are checked individually by our staff and CE conformed.

We take action to minimize the negative ecological impact of our industrial activity.

We are involved in energy saving policy.

**The maximum quality and service value...
all in one machine!!!**

**WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER**

max steam
LAB MANUFACTURER since 1966



for who believes and chooses reliability, safety, and easy use

The steamer with first rate quality spares and performance:

- Stainless steel boiler mm4
- Works with ordinary water
- Pressure switch
- Light level indicator
- Top professional electrovalve
- Steam flow control
- The most professional insulation system bandaging.

Energy saving policy

- Transformer
- Buzzer
- Cooper coated heating element with brass flange.
Energy saving policy
- Thermostat saver heating element
- The most original **MSgun** to distributes steam flow
- Top security with safety cap, fuse, pressure switch and transformer

We have our original gun and now with different way holding support.
The steamers are available with gun or nozzle and actived by button or pedal.
Manual and automatic way and now also with possibility to regulate and set the pressure level.

Eco line also with more many spares that respect the energy saving policy





MS 2

Steam cleaner ideal for a speedy use in little laboratories and dentists.

Platinum edition available also with **MSgun**.

Also available with box for ceramist works.

2lt 3 bar 1000W kg8 16x30x28cm

MS 3

The most used cleaner.

Steamer machine perfect for all the applications.

Unique in all the world !

Available also in Platinum Edition with 3lt 3 bar

1300W kg 12 18x30x24cm



MS 4

Steam generator with particular and symbolic design.

Highly reliable.

4lt 4 bar 1500W kg18 21x36x39cm

**WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER**

max steam
LAB MANUFACTURER since 1966

MS SMILE

Happy every time your sprinkling steam!
Stainless steel cover, versatile, safety, easy to use and service through innovative system of design.
Electrical components separated from the others for top security use and service
5lt 4,5 bar 1500W kg13 26x31x34cm



MS 5

The original steamer, the original since 1966.
Top machine and service performance
Perfect for medium and big laboratories for every day use of steam in a professional way.
5lt 4,5 bar 1500W kg14 22x31x37

MS LUXURY

Very precious steamer to give value to your job!!!

Reliable, experienced steamer
5lt 4,5 bar 1500W kg14 22x30x34cm





MS SPRINT

Fully attention to energy saving policies, many more spares to maintain heat inside and save energy.

For your Eco lab!!

Avalabe also in stainless steel version

5lt 4,5 bar 1500 W kg15

21x44x36cm



MS INFINITY

Just fill in always the water and work does not stop.

Infinite possibilities to use the clean generator

Never-ending steam!!

Gifted of electronic system control to set the water level.

Acoustic alarm when the water run out.

3lt + 1,5lt 4 bar 1500W kg13 26x31x34cm

MS UP&DOWN

Move and set the pressure for your right steam flow by electronic system!!!

You can go up & down with temperature system so means that you can change the pressure in every time.

Works with your right favourite steam flow

5lt 1500 W 5,5 bar 26x31



WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966

MS FULLSMILE –automatic

Connection with aqueduct.
Electronic feeler control to set the water level.
Ideal for big laboratories.
Professional use with steam by “MSgun”
and hot water by nozzle.
3lt 5bar 2000W kg17 26x31x34,5cm



MS POWERSMILE automatic

Connection with aqueduct. Electronic feeler control to set water level. You can also set the pressure right for your works and change it every time.

3 outlet: MSgun with steam, nozzle with water, **nozzle2** water and steam.

Powersteam in your lab!

3lt 6 bar 2000W kg 17
26x31x34,5cm



MS MasterSteam

Big steamer ,automatic version with main water connection. Electronic system , powerfull heating. Fitting box for clean with hot water and steam by jet erogator.

Stainless steel cover for professional and precious job.

2400W 5 bar

On request addition option air system





MS-MAXFLEXINJECTION


Injection flex press. The one made all in Italy.

With easy and professional electronic system you can set temperature and time for each kind of material and works. Available also a full kit with aluminium flask (single or double), duplication no-bump flask, micromotor, finishing material kit. Bio-compatible material with different kind of colors and transparent too. Real quality material bio compatible with oral mouth tissue

Material be relined. Rebasings units

Digital control, automatic injection.

Fast & easy process.

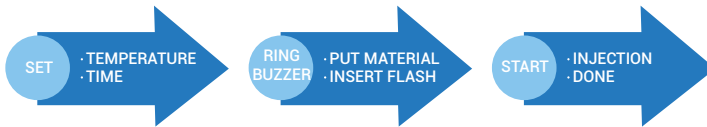

Select your data, you can program:

Automatic injection, simple easy, fast performance, save time and material.

Powerfull pression to obtain right precise, accurate and expanded works.

Micromotor 25000rpm reverse mode 80W


Time: 0-60 temperature: 0-300° or version 0-400°
6 bar

20x32x45cm kg 22 600W



WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966

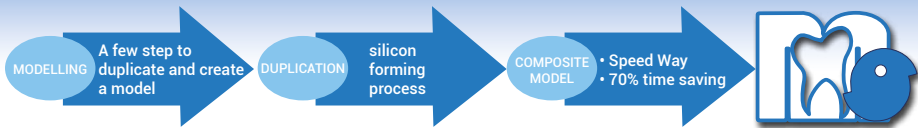

RESINE FLASK



ALUMINIUM

Single or double flask for injection device.
 Strong and precise.
 Aluminium material
 Available also in double version.

DUPLICATION

Stone duplication. Save time and material
 Many size. Anatomical band. Accurate and exact duplication
 Macro model to replicate every kind of dimension with extreme precision. No-bump material



DUPLICATION-COMPLETE DENTURE


Processing the denture using pouring technique.
 Bump free.
 Metal free.

👍 More production

👍 Saiving time

👍 Top performer

COMPRESSION MODELING TECHNIQUE for composite and acrylic resin

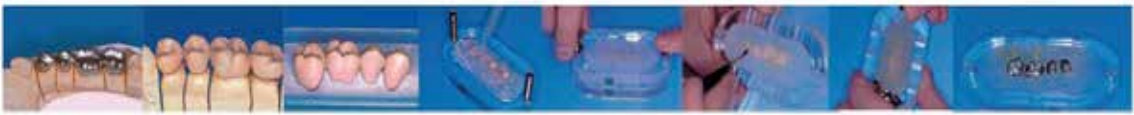



MS HMC

Heating, Melt, Cleaning... flask, cilinder, tool.
You can set time and temperature for many operations in the cleaning process
40x15x28cm W600 temp: 0-130°



MS CLEANING BOX

Wide range of cleaning box in stainless steel with removable floor. Plexiglass box with hand support and different dimensions
40x30x61cm 24x22x20cm 40x35x24cm 30x25x30cm

Now available also with LED LAMP!!!

Box support to cleaning turbine, with steam air and oil.



MS HY PRESS

Working pressure 220 bar, max pressure 400 bar , force 7500kg, 75000N , Max force 15000kg-150000N

MS BASKET

Make in aluminium,
strong use every day.



WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966



MS-MIXER

Last generation mixer, with low noise motor reducer and electronic system to selection program, time, speed and sense

Set programs in manual way or use 30 automatic mixing cycles (10 each plaster, coating and alginate)

Safety spares and buzzer stop

Available with wall or table support.

Wide range of cups with 200cc, 500 cc, 1000 cc and 2000cc

MS-VIBRATOR

Vibrator with serial electronic control intensity used for a wide range of drippings.

Uniform and constant steady vibration

Plate dimension 23x14,5cm kg 5 150W



MS-QUEEN

Pressure polymerisator for crowns and bridges. Safety with double stainless steel chambers 14cm.

Top precision microprocessor control polymerizator time (1-99min) and temperature (40-125°), adjustable for hot and cold system 1200W
35x35x26cm 13 kg



MS-3LIGHTS

Photopolymerization all kind of composite and resin. Works with 3 different light: 1 halogen and 2 leds about 45 degree. Also used flexible way (double led or without halogen)

Chamber 15x16x10cm 8kg 32x38x40 cm 150w
380-550Nm

MS-HANDY LIGHT

Compact and functional device to reduce processing polymerisation time, used with all composite coating applied before the final baking in the oven. Fixed the composite coat. Dichroic 100W activated by a practical aluminium button lamp 12 V 12x16x24,5 1,4 Kg 550-550nm



MS-BRIGHT

Compact and functional device to reduce processing polymerisation time, used with all composite coating applied before the final baking in the oven. Fixed the composite coat. Dichroic 100W activated by a practical aluminium button lamp 12 V 12x16x24,5 1,4 Kg 550-550nm

WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966



MS ROLLFURNACE

Cylinder furnace , easy ,compact ,functional and durable.It is possible set a final temperature to hold time at the end cycle to keep the cylinder in constant temperature.Electronic board control for program, time and temperature. Wide program range, each one contains 4 phases with 3 variables function (temperature, increase rate and holding time). Also the opportunity to set ignition delay and fumes aspirations. 3 different chamber size
1500W-2000W-2800W 1100° 28-38-45KG
CHAMBRE 16X16X10-18X23X12-23X30X15

MS WARME e MS WARME UP

New generation ceramic furnace for each kind of ceramics. Electronic system board to set program, update and memorize cycle in every phases works (100 program with 15 variables). It is possible set the vacuum level and regulate the cooling speed . 30kg 40x44x50 10A 1200W Available also with injection option (20 program) working by compression air. Used to fire all kind of ceramics, heated the solid ceramic up to its plastic status. Evolution furnace MS WARME UP
CHAMBRE 16X16X10-18X23X12-23X30X15





MS TRIMMER

Historical and simple functional trimmer.
Bakelite disk and water version
Simple and compact
1400-2800rpm 600W disc d14,7cm



MS -TRM 18

Trimmer dry and water version. With disk in bakelite or diamond. diameter 10" hole 1"
Powerfull engine, strong equipment.
Speed 14000rpm 1200W 24 kg
34x44x35cm



MaxSteam "Joiner" Saw

Extremely precise working thanks by innovative stumps you can cutting the models form the bottom.
Laser signal show the cut position and dual sensor ensures maximum ergonomics works in total safety way. 110V-220V 50/60Hz 200W 9000 RPM 26.5mm
Laser 2 10Kg 29x32x17



MaxSteam Saw

Saw to cut any kind of pattern with high precision. Saw with magnetic table and diamond disc.
Dimension 200x420x340 10KG Disc 80x0.2mm 60W
2800RPM



WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966



MS POLISHING

Polishing machine with double speed, equipped with two truncated cones.

Safety metal protection is available

1400-2800rpm 500W kg 7

MS SANDBLASTER

1 or 2 or 3 module machine

With pre -arranged power suction.

Pressure regulation

Available also model with RECYCLING and ECOLERSON

Movable nozzle support, halogen lamp 60W

10/25W 0.5-4 bar kg 13-16



MS WAXSHAPE

Double spatulas with electronic system for wax modeling. Each spatula can split On/off function and temperature control ; move in temperature range at 50°C to 200°C.

The electronic display show the temperature and it is possible move up &down degree by degree.

Easy and handy grips of hand pieces .Equipped with two separated hand pieces holders and with 5 interchangeable tips . 20W 1 kg

MS KSTONE

Natural stone for ceramics process work to maintain constant the humidity in natural and symbiosis way with labor room

All materials stay dry with two separated hand pieces holders and with 5 interchangeable tips.

20W 1 kg





MS FLASH

Professional laser pointer allows a rapid, precise and perpendicular drilling of plaster models.

High standards of safety of the powerful engine, activated by a micro-switch under the sliding work surface.

Adjustable drilling depth and with universal chuck. It is possible mount any type of bur and you can work all types of plaster, even those extra-hard, and plastics such as Plexiglas, PVC and resins.

90 W 5.000 rpm 15x14x31 cm 5 Kg

MS UFO

3 curing modes: intense, pulse and incremental light (20sec)

Every 10 sec Sound effect "beep"

Metal structure: very resistant 1300 mW/cm²

Long life lithium battery (180-195 application)

Dimension Length: mm 270x 24 Thickness



MaxSteam FINISHER

Finisher used for any type of resin, plaster or plastic material. Safety measures, flexible operation by different levels of control.

Dimension 210x280x290 350W 11KG

Bur connection d 6mm 28000 RPM

WE ARE EXPERIENCE
WE MAKE QUALITY
WE JOIN TOGETHER

max steam
LAB MANUFACTURER since 1966



MS CHEW

INNOVATIVE EQUIPMENTS

Thermoforming device to realize highly precise thermoformed performance .

Used by dentist and laboratories as it is very simple use and precise works.

Full automatic process thanks a professional board and with a silent and powerful vacuum pump & piston

- You can use round and square discs Ø 125 mm even thicker from 0.5 up to 5 mm
- Adjustable time selection
- Digital control panel-LCD display
- SMD technology
- Security system; Fault button
- The machine doesn't require compressed air or external air sources
- Removable basket for beads
- Well for realizing casing patterns
- 24x28x33 5,8kg 450W 250mbar
- ALL MADE IN ITALY. Check one by one



MS SCAN ART

MS ESTHETICS SMILE MATERIALS



MS CAD DISK



COLD LED DRYER



MICROMOTOR 40000rpm with LIGHT&SPRAY



CAD FINISHING TANKS

MS PHOTOPOLYMER MATERIAL

MS ACRYLIC

MS RESIN

MS SILICON GEL

MS SILICON

MS PHOTOSENSITIVE RESIN



100% ITALIAN ORIGIN

MODEL PREPARATION



MODELLING



DUPLICATION



FINISHING

